



Blood and fluid warming device
– safe and reliable

Barkey S-line

Barkey S-line is a one-touch-control infusion and blood warmer. The device helps to accelerate the healing process of patients and conserves hospital's budgets.

Benefits at a glance

- Warms infusion solutions, FFP and blood up to the patient
- Designed for any standard I.V. set
- One preset temperature for easy operation
- One-touch control
- No additional disposables required
- Displays actual temperature of the warming sleeve
- Use of biocompatible materials
- Visual and acoustic alarm in case of error
- Available in two different diameters (5 mm, 7 mm) and lengths (1.2 m, 1.5 m)
- Can be applied for the return of dialyses fluids


Barkey S-line

Safety

- Integrated temperature sensors
- Overtemperature alarm
- Smooth surfaces, therefore simple to wash and disinfect
- Operational safety thanks to preset target temperature value

Prevention of hypothermia

Warmth is the key to patients' wellbeing. Barkey's inexpensive infusion and blood warmer S-line helps to ease the treatment for all patients and accelerates the heating process. As a result care personnel has less work and treatment costs are reduced.

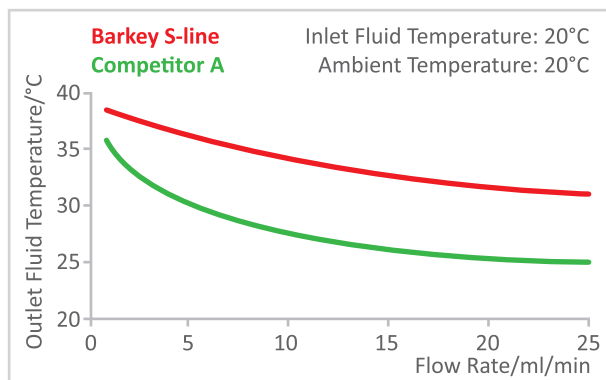

Product	Barkey S-line
Mains supply	100-240 V, 50/60 Hz
Power consumption	115-139 VA
Set temperature	+37°C oder +39°C
Length warming sleeve	1.2 m, 1.5 m
For I.V. set diameter	5 mm or 7 mm
Access line	3 m
Weight	approx. 1.3 kg

MDD classification	Ila
Protection class	I, protective earthing
Moisture protection	IP X2
Degree of protection	BF, defi-proof / internal and external treatment to the patient
Certification	CE 0123
Barkey is certified	DIN EN ISO 13485